



INSTRUCTIONS


TISSUE SAMPLER

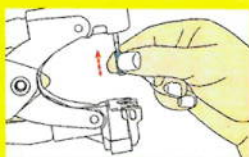

Load the blue cutter tip onto the tagger pin, holding it from the base, not from the tip.

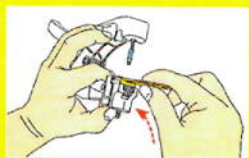

Push the metal plate up and slide the sampler tag underneath it. Clip the sampler flag onto the tagger lug. Check that the cutter tip is correctly aligned with the female by closing the tagger until the cutter tip reaches the female tag.

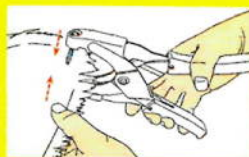


Move the tag away from the ear. Compress the tagger handles completely to take the sample.

TISSUE SAMPLER+TAG



Slide the male tag onto the tagger pin. Position the male flag on the left or right side.

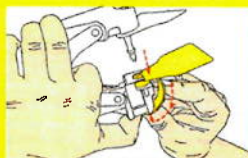

Slide the female tag into the tagger, underneath the metal plate. Place the container flag onto the tagger lug. Position the female flag to the side. Close the tagger and check the alignment between cutter tip and female tag.



Position the tag in the middle of the ear with the male tag facing the back of the ear. Drive the tagger all the way through the ear to feel a double click.

The Z Tags Tissue Tagger has been designed to apply Z Tags Tissue Sampler, and Tissue Sampler and Tag products, Z Tags two piece range and the Alpha Line range.

For more information on other quality Datamars brands and products, visit www.datamars.com or www.ztags.com

Manufactured in China.